



DVE 5500-D-GT



DVE 5500-D

The medical ceiling supply unit for a wide range of applications

- DVE 5500/5600-S - rigid horizontal unit
- DVE 5500/5600-D - rotating horizontal unit
- DVE 5500/5600-D-GT - with equipment carrier
- DVE 5500/5600-BSP-S - rigid vertical unit
- DVE 5500/5600-BSP-D - rotating vertical unit

modul
technik GmbH

*medical
supply systems
and equipment*

Ceiling Supply Unit - horizontal model

■ product description

- DVE 5500 (90°-, 70°- and 20°- level)
- DVE 5600 (90°- and 60°- level)

The ceiling supply units, DVE 5500/5600 is equipped on two/four levels with installation points for high and low voltage equipment and supply points for medical pressurized gases, vacuum, anaesthetic gas system (AGSS) and supply points for the provision and supply of air for the operation of surgical tools (air motor). The media are fed via the intermediary ceiling in to the rigid or rotating support columns.

The 90°-level is to be equipped on both sides with equipment rails **DIN EN 12218**.

The rotational movement of the DVE +/-165° is fixed by a stopper. The restriction of the rotational movement is to be adjusted using an adjustable friction brake.

The Ceiling Supply Units **DVE 5500/5600** are manufactured in compliance with recognized technological principles, and in particular, according to the current regulations for medical supply units in accordance with **DIN EN ISO 11197:2005**.

The available width for the receipt of the equipment consoles is 735 mm. The ceiling supply units are available in the following widths: **435, 635 and 735 mm**.



DVE 5500-D

fig.3



DVE 5600-D-GT
(with equipment carrier and console)

fig.4



DVE 5500-D-GT (with equipment carrier)

fig.5

Ceiling Supply Unit - vertical model

■ product description

- The basic support heads from the **DVE 5500-BSP-S/-D** and **DVE 5600-BSP-S/-D** range can be fitted with elements on three/five levels of installation.
- They are available in standard lengths of **600**, **1025** and **1325 mm**.
- The equipment consoles can be affixed when required to the front facade inside the guide rails, which are integrated into the basic support head, and can be adjusted to the exact required height.
- modul technik GmbH** has a wide range of accessories, which can be added as modules to ensure that the individual requirements of each work space are met. Why not ask for our accessories catalogue?
- Support pipe for infusion equipment**
The basic support heads of the DVE 5500- BSP / DVE 5600-BSP range are to be equipped at the 90°-level with stainless steel support pipes, Ø 30 mm, to receive swinging infusion pipes in standard lengths: **1000**, **1200**, **1400** and **1600 mm** and other equipment components.



DVE 5500-BSP-D with a console

fig.6



DVE 5500-BSP-D



fig.7





DVE 5500-BSP-D with consoles

fig.8

equipment level vertical

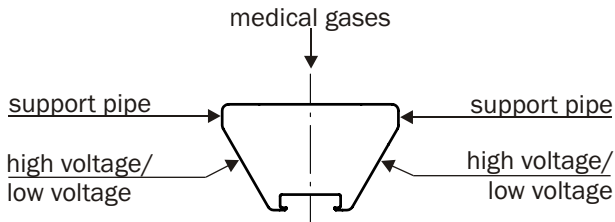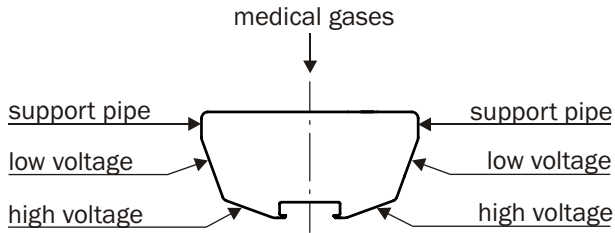


equipment level horizontal

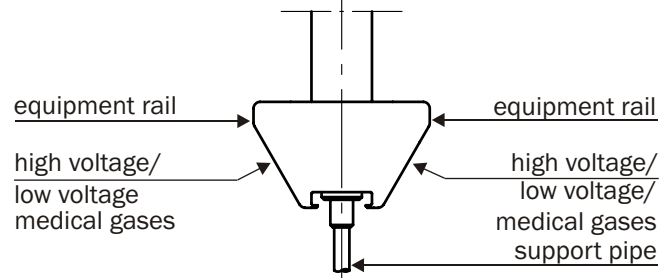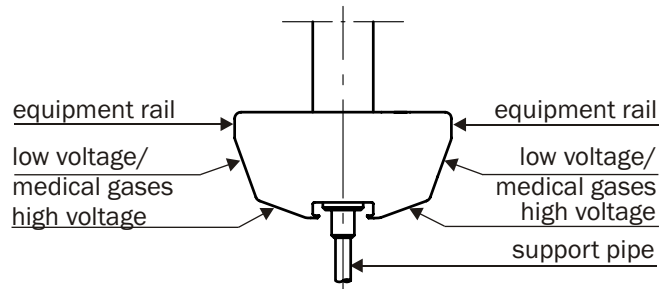

Our standard RAL-colours (other colours available on request)



RAL 9010
pure white

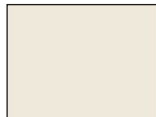


RAL 9001
cream



RAL 5014
pigeon blue



RAL 9002
grey white

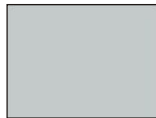

RAL 7035
light grey

(The colours featured in our brochure match as closely as possible the RAL chart)


The Quality Management System of **modul technik GmbH** is certified according to **DIN EN ISO 9001:2000**, **DIN EN ISO 13485:2003** and Appendix II of guideline 93/42/EWG (medical products).

modul
technik GmbH

*medical
supply systems
and equipment*

Rudolf-Diesel-Straße 5
D-56410 Montabaur
(Industriegebiet)
Tel.: (+49)2602/9449-0
Fax: (+49)2602/9449-11
e-mail: info@modul-technik.de
www.modul-technik.de

distributed by:

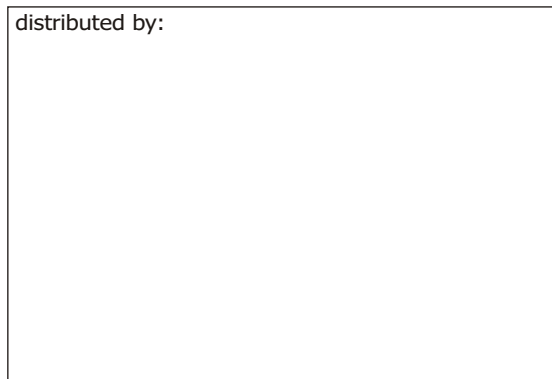

DVE 5500/5600-10/06-E2

CE 0366